See discussions, stats, and author profiles for this publication at: https://www.researchgate.net/publication/231205834

Determination of Picogram Quantities of Vanadium in Calcite and Seawater by Isotope Dilution Inductively Coupled Plasma Mass Spectrometry with Electrothermal Vaporization

ARTICLE in ANALYTICAL CHEMISTRY · JANUARY 1996		
Impact Factor: 5.64 · DOI: 10.1021/ac950859c		
CITATIONS	READS	
36	30	

2 AUTHORS, INCLUDING:



47 PUBLICATIONS 822 CITATIONS

SEE PROFILE

Determination of Picogram Quantities of Vanadium in Calcite and Seawater by Isotope Dilution Inductively Coupled Plasma Mass Spectrometry with Electrothermal Vaporization

David W. Hastings* and Steven R. Emerson

School of Oceanography, WB-10, University of Washington, Seattle, Washington 98195

Bruce K. Nelson

Department of Geological Sciences, AJ-20, University of Washington, Seattle, Washington 98195

We developed a method to measure picogram quantities of vanadium in calcite and seawater by isotope dilution (ID) inductively coupled plasma mass spectrometry using electrothermal vaporization (ETV) to introduce the sample into the plasma. \overline{A} 50V isotope spike enriched to 44 atom % was equilibrated with samples, followed by chemical purification by cation exchange chromatography. Samples were introduced into the ETV unit with a Pd modifier and heated to 1000 °C. This quantitatively eliminates the ClO⁺ isobaric interference with V at m/z 51 for solutions up to 0.5 N HCl. The procedural blank was 0.27 pg of V. Corrections for ⁵⁰Ti and ⁵⁰Cr, which interfere with the ⁵⁰V signal, were made by measurement of ⁴⁹Ti and ⁵³Cr. These isobaric interferences and variable ArC levels were the limiting sources of error in the ID measurement and diminished the detection limit to 6 pg of V. The detection limit for nonisotope dilution applications was 0.3 pg of V. Measurement precision on the same sample of dissolved calcite over the course of one run was $\pm 3\%$ (1 σ). Accuracy was confirmed by determination of V standards in CaCO₃ and by comparative measurement with ID thermal ionization mass spectrometry and graphite furnace atomic absorption spectroscopy.

Calcitic shells of foraminifera, unicellular marine animals that grow in surface waters and in ocean sediments, have concentrations of certain elements that are determined by the seawater concentration of the water in which the shell was formed. The determination of trace elements such as Cd, Ba, and Li in foraminiferal calcite has proven to be a powerful tool to deduce the history of ocean circulation patterns, nutrient concentrations, and oceanic inputs.^{1–4} The development of new foraminiferal tracers including B, F, U, Mg, and V is an active area of oceanographic research which has shed new light on paleoclimatology and variations in the global carbon cycle.^{5–8}

Recent interest in the geochemistry of vanadium has arisen from its unique redox behavior 9,10 and the observation of enriched concentrations in anoxic sediments, black shale, and crude oil. 11,12 Because V deposition is redox sensitive, its concentration in seawater is most likely sensitive to the oxidation state of sediments. The calcitic shells of foraminifera accurately record seawater concentrations of $V,^{13}$ so the V content of the fossil shells of these microscopic oceanic organisms potentially can serve as an indicator of changes in past ocean anoxia. 14

A typical 5 mg hand-picked sample of these microscopic shells yields only 50-200 pg of vanadium. Existing analytical methods for the determination of vanadium at such low levels are inadequate because of poor sensitivity, calcium matrix effects, and isobaric interferences. Table 1 summarizes methods to determine the concentration of V in seawater and other low-level environmental samples. The most common methods use either instrumental neutron activation analysis (INAA) or graphite furnace atomic absorption spectrometry (GFAAS). A preconcentration step is used, which typically involves coprecipitation with ferric hydroxide,15 cobalt APDC,16 or APDC and activated carbon.17 Addition of coprecipitants necessarily increases the procedural blank. Potential variability in the recovery for the preconcentration results in inaccurate determinations for several of these methods. As a result, significant variations have been obtained for surface seawater samples. 18,19

^{*} Currently at Department of Oceanography, University of British Columbia, 6270 University Blvd., Vancouver, BC V6T 1Z4, Canada.

⁽¹⁾ Boyle, E. A. Annu. Rev. Earth Planet. Sci. 1992, 20, 245-287.

⁽²⁾ Delaney, M. L.; Bé, A. W. H.; Boyle, E. A. *Geochim. Cosmochim. Acta* 1985,

^{45, 326 336.} (3) Lea, D. W.; Boyle, E. A. *Nature* **1990**, 347, 269–272.

⁽⁴⁾ van Geen, A.; Luoma, S. N.; Fuller, C. C.; Anima, R.; Clifton, H. E.; Trumbore, S. Nature 1992, 358, 54–56.

⁽⁵⁾ Sanyal, A.; Hemming, N. G.; Hanson, G. N.; Broecker, W. S. Nature 1995, 373, 234–236

⁽⁶⁾ Opdyke, B. N.; Walter, L. M.; Huston, T. J. Geology 1993, 21, 169-172.

⁽⁷⁾ Rosenthal, Y.; Boyle, E. A. Geochim. Cosmochim. Acta 1993, 57, 335-346.

⁽⁸⁾ Russell, A. D.; Emerson, S.; Nelson, B. K.; Erez, J.; Lea, D. W. Geochim. Cosmochim. Acta 1994, 58, 671–681.

⁽⁹⁾ Calvert, S. E.; Pedersen, T. F. Mar. Geol. 1993, 113, 67-88.

⁽¹⁰⁾ Wehrli, B.; Stumm, W. Geochim. Cosmochim. Acta 1989, 53, 69-77.

⁽¹¹⁾ Breit, G. N.; Wanty, R. B. Chem. Geol. 1991, 91, 83-97.

⁽¹²⁾ Lewan, M. D.; Maynard, J. B. Geochim. Cosmochim. Acta 1982, 46, 2547– 2560

⁽¹³⁾ Hastings, D. W. Vanadium in the Ocean: A Marine Mass Balance and Paleoseawater Record. Ph.D. Thesis, University of Washington, Seattle, WA, 1904

⁽¹⁴⁾ Emerson, S. R.; Huested, S. S. Mar. Chem. 1991, 34, 177-196.

⁽¹⁵⁾ Weisel, C. P.; Duce, R. A.; Fasching, J. L. Anal. Chem. 1984, 56, 1050– 1052.

⁽¹⁶⁾ Shiller, A. M.; Boyle, E. A. Earth Planet. Sci. Lett. 1987, 86, 214-224.

⁽¹⁷⁾ Middelburg, J. J.; Hoede, D.; van der Sloot, H. A.; van der Weijden, C. H.; Wijkstra, J. Geochim. Cosmochim. Acta 1988, 52, 2871–2878.

⁽¹⁸⁾ Collier, R. W. Nature 1984, 309, 441-444.

Table 1. Comparison of Methods for the Determination of Vanadium in Seawater and Other Environmental Samples

	1 5773
Co-APDC GFAAS 16 500 pg 5 80 seawa ID-TIMS 36 90 pg 12 biolog ID-TIMS a 18 pg 5 93 calcite ID-ICPMS-ETV a 6 pg 3% 93 calcite	ted [V] _{sw} is 30% low ater analysis gical samples e, seawater analysis e, seawater analysis determination

Determination of vanadium by inductively coupled plasma mass spectrometry (ICPMS) is complicated by the presence of ³⁵Cl¹⁶O⁺ in the plasma. In samples with even small amounts of chlorine (0.01%), these polyatomic ions overlap the spectra of the elements of interest. A wide variety of methods have been developed to correct or reduce the ClO interference. One of the simplest is to determine the 35ClO concentration by using the known isotope ratio of 35Cl/37Cl and monitoring the 37ClO peak at m/z 53.²⁰ This is feasible only if Cr levels in the sample are low and if relatively low levels of Cl (<1%) are present. Alternative sample preparation methods which eliminate the use of HCl^{21,22} are a first-order solution to reducing some polyatomic interferences. On-line separation methods to remove Cl species which cause spectral interferences have been successfully employed.²³ Cryogenic desolvation has been found to be effective in reducing the levels of interfering oxides and chlorides by removing the water and HCl in the sample.²⁴ Anodic stripping voltammetry interfaced with ICPMS eliminates the ClO+ interference, but only at levels up to 0.1% chloride.25

A variety of mixed gas plasmas has been used to attenuate polyatomic interferences; for example, N_2 in the outer (coolant or plasma) gas²⁶,²⁷ and N_2 , Xe, O_2 , or CH₄, mixed with Ar in the aerosol (nebulizer) gas.^{28–30} High-resolution mass spectrometry to resolve the spectral interferences is an effective yet costly technique,³¹ since a resolution of 30 000 is necessary to resolve ClO⁺ at m/z 50.964 and V⁺ at m/z = 50.944. The wide variety of approaches used in resolving or attenuating interferences in ICPMS are included in a comprehensive review.³²

In this paper, we describe a method to determine vanadium in calcite and in seawater by isotope dilution (ID) ICPMS using electrothermal vaporization to eliminate the volatile chlorine species. In addition to removing the isobaric polyatomic ions, ETV dramatically increases the efficiency of analyte transport into the plasma compared to conventional nebulization, from <2% to

60-80% This improves the detection limit by at least an order of magnitude. Use of isotope dilution and no chemical preconcentration results in highly reliable determinations and low procedural blank.

EXPERIMENTAL SECTION

Standards and Solutions. Solutions were prepared using distilled, deionized (DDI) water. Reagent grade HNO3 and HCl were doubly distilled in a quartz still. For ultratrace applications, acids and water were purified further by triple and quadruple distillation in Teflon (FEP) sub-boiling stills and stored in acidleached Teflon bottles. Ultrapure H₂O₂ (Ultrex, J. T. Baker) was stored in the dark at 4 °C. Ten parts-per-million Pd solution was prepared by dissolving 99.9% Pd(NO₃)₂ (Aesar/Johnson Matthey, Ward Hill, MA) in sub-boiled HNO3. Teflon (PFA) beakers (Savillex Corp., Minnetonka, MN) were boiled in 8 N HNO₃, rinsed in DDI H₂O, and refluxed for 24 h on a heated clean table with 1 N doubly distilled HCl, followed by refluxing for 24 h with $2\times$ sub-boiled water. All plasticware was leached in 1 N HCl at 50 °C for at least 12 h; polypropylene autosampler cups were immersed in hot 1 N HNO3 for a minimum of 3 days and then rinsed with DDI water and dried immediately before use.

Enriched vanadium-50 isotope spike was obtained as V_2O_5 from the National Institute of Standards and Technology, originally prepared at Oak Ridge National Laboratories (ORNL). Approximately 0.5 mg of V_2O_5 was dissolved in ultrapure HNO $_3$ and diluted with 0.1 N HNO $_3$ to meet the requirements of the measurement. Original isotopic analysis performed at ORNL was stated as 44.1% ^{50}V (±0.1), 55.9% ^{51}V (±0.1), or a ratio of $^{50}V/^{51}V$ = 0.789. The isotope ratio of the spike was determined in this laboratory by thermal ionization mass spectrometry (TIMS) and by ETV-ICPMS. The V concentration of the spike solution was determined by isotope dilution TIMS by comparison to a gravimetrically prepared V single-element ICP standard accurate to within \pm 0.3%.

Vanadium standard for vanadium analysis in calcite was prepared by cleaning 5 g of bulk foraminiferal calcite, dissolving this in concentrated HNO3, and diluting to $\sim\!150$ mL with DDI H_2O . This standard was analyzed to determine the relative accuracy and precision of both the TIMS and ICPMS analyses; the standard was analyzed every nine samples. Replicate measurements of the standard were analyzed to determine the relative precision of the measurement. A CaCO3 standard prepared from ultrapure CaCO3 (99.9965% Puratronic grade, Aesar/Johnson Matthey) and spiked with a known amount of the V single-element ICP standard was also used to determine the accuracy of each method.

⁽¹⁹⁾ Jeandel, C.; Caisso, M.; Minster, J. F. Mar. Chem. 1987, 21, 51-74.

⁽²⁰⁾ Henshaw, J. M.; Heithmar, E. M.; Hinners, T. A. Anal. Chem. 1989, 61, 335–342.

⁽²¹⁾ Jarvis, I.; Jarvis, K. E. Chem. Geol. 1992, 95, 1-33.

⁽²²⁾ McLaren, J. W.; Beauchemin, D.; Berman, S. S. Anal. Chem. 1987, 59, 610-613

⁽²³⁾ Goossens, J.; Dams, R. J. Anal. At. Spectrom. 1992, 7, 1167-1171.

⁽²⁴⁾ Alves, L. C.; Allen, L. A.; Houk, R. S. Anal. Chem. 1993, 65, 2468-2471.

⁽²⁵⁾ Pretty, J. R.; Blubaugh, E. A.; Caruso, J. A.; Davidson, T. M. Anal. Chem. 1994, 66, 1540-1547.

⁽²⁶⁾ Craig, J. M.; Beauchemin, D. J. Anal. At. Spectrom. 1992, 7, 938–942.
(27) Wang, J.; Evans, E. H.; Caruso, J. A. J. Anal. At. Spectrom. 1992, 7, 929–

⁽²⁷⁾ Wang, J.; Evans, E. H.; Caruso, J. A. J. Anal. At. Spectrom. 1992, 7, 929-936.

⁽²⁸⁾ Smith, F. G. Anal. Chem. 1991, 63, 1458-1462.

⁽²⁹⁾ Hill, S. J.; Ford, M. J.; Ebdon, L. J. Anal. At. Spectrom. 1992, 7, 719-725.

⁽³⁰⁾ Lam. J. W. H.: Horlick. G. Spectrochim. Acta **1990**. 45B. 1313–1325.

⁽³¹⁾ Morita, M.; Ito, H.; Uehiro, T.; Atouka, K. Anal. Sci. 1989, 5,

⁽³²⁾ Evans, E. H.; Giglio, J. J. J. Anal. At. Spectrom. 1993, 8, 1-18.

Sample Preparation. The overall chemical preparation scheme involves picking single species of foraminifera from a ocean sediment core, removing contaminant phases via chemical and physical cleaning methods, and dissolving the cleaned foraminiferal calcite. This solution was then spiked with the enriched isotope tracer, and major cations (Ca²⁺ or Na⁺) and some isobaric interferences (Cr and Ti) were removed by cation exchange chromatography. The sample was then dried and taken up in a small volume of dilute nitric acid for injection into the ETV-ICPMS.

Cleaning Procedure. Foraminifera were obtained by washing and sieving bulk ocean sediment several times to remove clays and fines. Single species of foraminifera were hand-picked from the washed, dried, size-fractionated sample. Several V-rich contaminant phases associated with foraminifera must be removed from the foraminiferal calcite prior to analysis. Foraminifera were gently crushed to open the individual chambers and rigorously cleaned using a method modified from Boyle³⁴ to eliminate these phases. Following the cleaning procedure, all sample handling was performed in a HEPA laminar flow hood in a class 100 clean room. For a complete description of the cleaning procedure, see ref 13.

Cleaned samples were transferred to 6 mL Teflon beakers, subjected to a final 0.005 N HNO $_3$ rinse to remove remaining adsorbed metals, and dissolved in 600 μ L of 6 N HNO $_3$. Following dissolution, 500 μ L of each sample was transferred to another set of 6 mL Teflon beakers for V analysis, gravimetrically spiked with the enriched 50 V isotope tracer, and taken to dryness to ensure isotopic equilibration. The remainder of the sample was used for Ca, Mg, and Mn analyses. Dried samples were redissolved in 500 μ L of 0.01 N HCl and 0.3% H $_2$ O $_2$ and loaded onto prepared ion exchange chromatography columns.

Ion Exchange Chromatography. The cation exchange resin (AG50W-X8, 100-200 mesh, H^+ form; Bio-Rad Laboratories, Hercules, CA) was batch cleaned prior to use by soaking 150 g of resin in 750 mL of 6 N HCl for at least 1 h. This treatment was repeated six times and followed by repeated rinses with DDI water until pH ≥ 4 and stored in DDI water for at least 2 weeks prior to use. For 3-10 mg samples, 1.5 mL of hydrated resin slurry was poured into precleaned 10 mL polypropylene columns (Bio-Rad). The resin was cleaned in situ with one column volume of 1 N HCl and 0.6% H_2O_2 , followed by three column volumes of 0.01 N HCl and 0.3% H_2O_2 . The final column conditioning rinse was prepared with the cleanest acids and water available.

Cleaned, dissolved, and spiked samples were loaded onto the resin, eluted with 4 mL of 0.01 N HCl and 0.3% $\rm H_2O_2$, and collected in 6 mL Teflon beakers. A vanadium(V) peroxide complex forms which was eluted quantitatively while $\rm Ca^{2+}$ was retained on the column. Any reduced vanadyl should be oxidized to vanadate by the peroxide and eluted from the column. High titanium and chromium levels in the final sample will severely compromise both the accuracy and the precision of the vanadium measurement by isotope dilution. Both have isotopes at m/z 50 which are isobaric with $^{50}\rm V$, so it is essential to eliminate or reduce their concentration in the final sample. Titanium forms a complex with $\rm H_2O_2$ that is strongly retained by the cation exchange resin at the low acid strengths used; final $\rm ^{50}\rm Ti$ levels were typically $\rm ^{<2}$ pg. $\rm Cr(III)$ is

preferentially eluted at higher acid strengths (3 M HCl), but Cr-(VI) is eluted even at 0.01 N HCl; 50 Cr was usually $^{<6}$ pg in the final sample. Following the ion exchange step, samples were taken to dryness and redissolved in 0.16 N HNO₃. Processing for analysis by ICPMS and TIMS was identical except that sample size was substantially larger for the TIMS measurement (30–40 mg compared to 3–10 mg of calcite), so the resin volume was increased to 4 mL.

Chemical yields for the ion exchange step were determined by spiking the standard foraminifera mixture with the enriched isotope spike following the column chromatography. The calculated value was compared with the value determined by spiking the sample beforehand. Recovery of V from the ion exchange step varied from 89 to 95%, with an average recovery of 93%. It should be noted that, assuming isotopic equilibration between spike and sample before any sample processing, variability in recoveries will not affect the final result.

Seawater samples (0.1–1 g) were processed in the same way as calcitic samples. After acidification to pH 2, they were spiked with the enriched ^{50}V isotope tracer, taken to dryness and dissolved in 500 μL of 0.01 N HCl and 0.3% H_2O_2 . Elution through the cation exchange resin eliminates the matrix interferences from the high sodium concentration and removes some of the trace Ti and Cr.

Thermal Ionization Mass Spectrometry. Vanadium was analyzed on a VG Sector thermal ionization mass spectrometer equipped with a Daly detector for signal amplification. Samples were dried down to a $2-3 \mu L$ drop, loaded onto a zone-refined, outgassed Re filament with a colloidal graphite suspension, and gradually heated to 1400-1450 °C.36 To reduce contamination, the graphite was leached twice in 6 N HCl, rinsed with DDI H₂O, and resuspended in $4 \times$ sub-boiled water prior to use. The signals at m/z 49.5, 49, and 52 were monitored in order to correct for background and isobaric interferences. 50Ti and 50Cr were estimated by measuring 49Ti and 52Cr and using the natural isotopic abundance ratios $^{50}\text{Ti}/^{49}\text{Ti} = 0.954$ and $^{50}\text{Cr}/^{52}\text{Cr} =$ 0.051 86 to determine the appropriate correction factor.³⁷ Two blocks of data with 10 individual ratio measurements of 5 s integration were taken for each sample. Total time to run each sample was \sim 50 min. Minimum sample size was 150 pg of V, or about 30 mg of calcite before cleaning. Typical observed ion currents were $1 \times 10^{-15} - 4 \times 10^{-15}$ A on 51 V.

Inductively Coupled Plasma Mass Spectrometry. A Perkin Elmer Sciex ELAN 5000 ICPMS with a Perkin Elmer HGA-600MS electrothermal vaporization (ETV) unit and an AS-60 autosampler was used. The conventional detector (channel electron multiplier) was replaced with a discrete dynode electron multiplier (ETP Scientific, Auburn, MA). The ICP operation conditions and ETV temperature program are given in Table 2. Ion optics, nebulizer flow, and torch position were optimized for maximum sensitivity and minimum noise for vanadium m/z 51 while the conventional pneumatic nebulizer was in place. Doubly charged ions and oxide levels, as represented by Ba²⁺/Ba⁺ and CeO/Ce, were monitored in the optimization process and kept below 2 and 3%, respectively. No readjustment of operating conditions was necessary after attaching the ETV unit. Signals were collected in peak height mode. Pyrolytically coated graphite tubes were used. To remove

⁽³⁴⁾ Boyle, E. A.; Keigwin, L. D. Earth Planet. Sci. Lett. 1985/86, 76, 135-

⁽³⁵⁾ Fritz, J. S.; Abbink, J. E. Anal. Chem. 1962, 34, 1080-1082.

⁽³⁶⁾ Fassett, J. D.; Kingston, H. M. Anal. Chem. 1985, 57, 2474–2478.

⁽³⁷⁾ Fassett, J. D.; Beary, E. S.; Xiong, X.; Moore, L. J. Anal. Chem. 1994, 66, 1027–1031.

Table 2. Operating Conditions for ELAN 5000 ICPMS-ETV

Inductively Coupled Plasma					
Ar outer flow	15 L min ⁻¹				
Ar nebulizer flow	$1.015~{ m L}~{ m min}^{-1}$				
Ar auxiliary flow	0.85 L min ⁻¹				
forward power	1.1 kW				
detector voltage	3.8 kV				
Mass	s Spectrometer				
lenses	tune to V at m/z 51				
sampling position	9 mm from load coil				
sampler cone	Pt, 1.1 mm diameter orifice				
skimmer cone	Pt, 0.9 mm diameter orifice				
dwell time	15 ms				
measurement time	675 ms				
resolution mode	normal (0.7 amu at 10% peak height)				
scanning mode	peak hop transient, 1 pt/peak				
signal profile	maximum (peak height)				
signal profile smoothing	11 point running average				
sample volume	$30-35 \mu L$ (in duplicate)				
modifier	200 ng of Pd as PdNO ₃				
	· ·				

Temperature Cycle for ETV

			argon flow (mL/mir	
temp (°C)	ramp (s)	hold (s)	internal	carrier
100	2	10	300	0
120	10	5	300	0
1000	5	15	300	0
2400	0	5	0	1000
20	1	5	0	1000
2650	1	5	0	1000
20	1	10	0	1000
	100 120 1000 2400 20 2650	100 2 120 10 1000 5 2400 0 20 1 2650 1	120 10 5 1000 5 15 2400 0 5 20 1 5 2650 1 5	temp (°C) ramp (s) hold (s) internal 100 2 10 300 120 10 5 300 1000 5 15 300 2400 0 5 0 20 1 5 0 2650 1 5 0

^a Drying steps 1 and 2 are carried out following addition of both Pd modifier and sample.

residual vanadium before sample analysis, graphite tubes were conditioned by heating to 2650 °C for 10 s repeatedly until tube blanks did not decrease significantly and were below 700 counts/s at m/z 51. The sample passes from the graphite furnace to the torch via a 63 mm long, 6.4 mm i.d. Teflon tube.

Palladium Nitrate. Addition of a carrier is essential to assure that the analyte does not condense on surfaces following vaporization in the graphite furnace. We used 20 µL of a 10 ppm PdNO₃ solution. By acting as a physical carrier of the vaporized analyte, the Pd assured adequate transport efficiency from the graphite tube into the plasma.^{38,39} Two hundred nanograms of Pd was determined to be the optimal amount for vanadium determination. The intensity decreases with lesser amounts, while additions larger than 200 ng of Pd do not significantly improve the signal. We explored several other modifiers including dilute seawater, lanthanum oxide, and the CaCO₃ present in the original sample. Ten microliter aliquots of seawater diluted 1:1000, which contains \sim 350 ng of Na and is intrinsically low in transition metals, was used. 40,41 While the seawater modifier increased the intensity relative to no added carrier, 200 ng of PdNO3 gave 30% more counts in peak height mode. Signal intensities for V standards using a La modifier were also inferior to those for Pd. The CaCO₃ in the original sample, typically 2 mg, was not a useful modifier; calcium suppressed the signal by 90%.

ETV Temperature Program. The Pd modifier was introduced into the graphite tube and heated gradually to 120 °C. After the tube was cooled briefly, 35 μ L of sample was introduced and slowly dried. The sample was charred at 1000 °C for 20 s to eliminate Cl and vaporized at 2400 °C to introduce it into the plasma. Maximum power heating was used in this step, with a heating rate of ~ 2000 °C s⁻¹. An optimal vaporization temperature of 2400 °C was chosen on the basis of obtaining the highest signal while prolonging tube life. At 2600 °C, the vanadium signal was 15% greater, but the tube life was significantly shortened. Three seconds before and during the vaporization step, the sample introduction hole in the graphite tube was closed with a graphite rod to prevent sample loss. During this period, the argon flow passes through the graphite tube, down the Teflon tubing and into the ICP torch, instead of being directed from the ends of the tube out of the sample introduction hole, as it was during the preceding dry and pyrolysis steps. Vanadium is a refractory metal and forms refractory carbides in the graphite furnace, 42 resulting in peak tailing and carryover of V from one sample to the next. To reduce this effect, a high-temperature cleaning step (2650 °C) was used following the vaporization step.

Intensity at m/z 49, 50, 51, 52, and 53 was measured with 15 ms dwell times in peak hop transient mode at normal resolution (Table 2). Background levels were determined at the beginning of the run by monitoring those masses for the Pd modifier only. Each sample was analyzed in duplicate. Total time for each duplicate analysis, including addition of Pd modifier and sample, drying, charring, vaporization, and cleaning, was <6 min.

Isotope Dilution Calculations. Isotope dilution allows highly accurate and precise determinations in elemental determination. Any variable sample loss following spike equilibration (e.g., during ion exchange chromatography, sample processing, or in the ETV) will not affect the final value. The concentration of vanadium in the sample spiked with the enriched 50 V tracer was determined from the isotopic ratio, $R_{\rm m}$, measured by mass spectrometry and calculated according to the following equation:

$$C_{x} = \frac{C_{s}W_{s}}{W_{x}} \frac{^{50}V_{s} - R_{m}^{51}V_{s}}{R_{m}^{51}V_{x} - ^{50}V_{x}}$$

where $C_{\rm x}$ and $C_{\rm s}$ are the concentrations of V in the sample and spike, respectively; $W_{\rm x}$ and $W_{\rm s}$ are the weights of the sample and spike; $^{50}V_{\rm x}$ and $^{50}V_{\rm s}$ are the atom fractions of the enriched isotope 50 V in the sample and spike, respectively; $^{51}V_{\rm x}$ and $^{51}V_{\rm s}$ are the atom fractions of 51 V in the sample and spike, respectively; and $R_{\rm m}$ is the measured 50 V / 51 V ratio of the spiked sample. $^{50}V_{\rm x}=0.0024$ and $^{51}V_{\rm x}=0.9976$ are taken from the natural abundances of the isotopes. $^{50}V_{\rm s}=0.444$ and $^{51}V_{\rm s}=0.556$ are taken from repeated measurements of the spike by TIMS. These values agree to within $\pm 0.3\%$ of the original isotopic analysis performed at OPNI

In isotope dilution calculations, the precision of the analysis depends to some extent on the value of the measured ratio, $R_{\rm m}$. An optimal ratio equal to the geometric mean of the spike and sample ratios, $R_{\rm opt}=0.002$, will theoretically minimize the error introduced by the propagation of errors in the above isotope dilution equation.⁴³ However, the uncertainty decreases with

(42) Wendl, W.; Müller-Vogt, G. Spectrochim. Acta 1985, 40B, 527-531.

⁽³⁸⁾ Gregoire, D. C.; Sturgeon, R. E. Spectrochim. Acta 1993, 48B, 1347-1364.

⁽³⁹⁾ Ediger, R. D.; Beres, S. A. Spectrochim. Acta 1992, 47B, 907–922.

⁽⁴⁰⁾ Hughes, D.; Chakrabarti, C. L.; Goltz, D. M.; Gregoire, D. C.; Sturgeon, R. E.; Byrne, J. P. Spectrochim. Acta 1995, 50B, 425–440.

⁽⁴¹⁾ Sturgeon, R. E.; Willie, S. N.; Zheng, J.; Kudo, A.; Gregoire, D. C. J. Anal. At. Spectrom. 1993, 8, 1053-1058.

increasing total counts for both isotopes, which results in a greater optimal ratio.⁴⁴ A target spiking ratio of $R_{\text{opt}} = 0.1 - 0.3$ was chosen.

Mass Bias. Error introduced by instrumental and matrixinduced mass discrimination is an important consideration in mass spectrometry and must be evaluated before these methods can be relied on to produce reliable and accurate results. Unfortunately, there is no good isotopic standard for vanadium. Since the abundance of ⁵⁰V is extremely low in natural samples (0.24%), it is not possible to accurately measure the natural isotopic ratio by TIMS or ICPMS. Therefore, mass bias effects must be estimated. In TIMS, the greatest source of mass fractionation is Raleigh distillation as the sample is ionized from the heated filament. Results from K and Rb in our laboratory indicate that mass fractionation was rarely >0.1% amu⁻¹ in the relevant mass range. The isotope ratio of the spike, as measured by TIMS in our laboratory and corrected for isobaric interferences 50Cr and ⁵⁰Ti, was 0.799 (± 0.001 ; 1σ). This value did not change within the measurement error over a period of 4 years. It agrees with the ratio measured at ORNL to within $\pm 0.3\%$, suggesting that mass fractionation by TIMS was consistently reproducible.

While such fractionation at the ion source does not occur in plasma mass spectrometry, significant mass bias effects are possible in the lens stack and within the quadrapole mass spectrometer. Instrumental mass discrimination has been measured at 1-2% amu⁻¹ for Os isotopes⁴⁵ and as large as 5-7% for ¹⁰B/¹¹B ratios. ^{46,47} Mass fractionation in the ELAN 5000 was estimated by comparing the isotopic composition of the vanadium spike as determined by TIMS with those measured by ICPMS. The enriched isotope spike was analyzed by ICPMS frequently over the course of several runs to assess the relative degree of mass bias within each run and between runs. The ⁵⁰V/⁵¹V ratio of this standard did not vary by more than $\pm 0.5\%$ amu⁻¹ during any single run, indicating that machine drift was not a significant concern. However, the isotopic ratio of the spike determined by ETV-ICPMS was up to 5% lower than the value determined by TIMS. A mass fractionation correction factor, α , is included in the calculation to determine the vanadium content and provide consistency with TIMS analysis: $\alpha = (R_{\rm m}/R_{\rm spike} - 1)/\Delta m$, where $R_{\rm m}$ is the measured ratio, $R_{\rm spike}$ is the ratio determined by TIMS, and $\Delta m = 1$ is the mass difference between the two isotopes ⁵⁰V and 51V. Given a mass bias of 3%, the final calculated concentration for a 200 pg V sample using a 0.1 g spike was 5% higher using this correction factor. An alternative method to treat mass bias concerns would be to make the isotope spike calibration measurements by ID-ETV-ICPMS instead of by ID-TIMS, and use this value in the calculations. Since, to first order, spike calibrations should exhibit the same mass discrimination as samples, any bias is canceled in concentration measurements, and no correction would need to be made.

RESULTS AND DISCUSSION

Blanks. Low-level vanadium analyses are possible only if reagent, water, and procedural blanks are consistently at extremely

Table 3. Vanadium Content of Water, Reagent, and Process Blanks

reagents	V concn (pg/mL)	amount used ^a (mL)	V total (pg)
H_2O	0.05	4.0	0.20
H_2O_2	5.62	0.01	0.06
HCl	0.13	0.04	0.01
total known so measured prod	0.27 0.27		
unknown sour	<dl< td=""></dl<>		

^a Four milliliters of 0.01 N HCl/0.3% H₂O₂ used for each sample. ^b Resin, containers, airborne contamination, etc.

Table 4. Isotope Abundances of Elements and Interfering Polyatomic Ions in the Relevant Mass Range for Vanadium Determination by Isotope Dilution

	m/z							
	47	48	49	50	51	52	53	54
V Ti Cr Fe	7.28	73.94	5.51		99.76	83.76	9.55	2.38 5.82
	PO^+	$^{48}\mathrm{Ca}^{+}$			ClO^+	$Ar^{12}C^{+}$	$Ar^{13}C^{+}$	

low levels. To avoid contamination, sample processing was kept at an absolute minimum, and reagent volumes were as low as possible without sacrificing sample recovery.

Blank levels for each of the reagents and the water were quantified by comparing the vanadium intensity at m/z 51 with those of known standards. Sample volumes of 20-30 mL were used in order to obtain a true blank value significantly above the limit of detection (DL), even though less volume was used in practice. The V content of one, two, and four combined procedural blanks was determined to obtain an accurate value for the blank. A regression of the V content versus the number of procedural blanks results in a slope of 0.27 pg of V per blank ($r^2 = 0.97$). Variability of this blank value was primarily dependent on the quality of the water used, although contaminated batches of resin occasionally contributed to a high blank. Abnormally high levels of chromium were also found in some resin batches. In these cases of higher V or Cr, the resin was recleaned by soaking in HCl, as described earlier. The total known sources of blank from the reagents and the measured procedural blank were equal (Table 3), which indicates that unknown sources of blank including V from the resin, containers, and environmental contamination were negligible.

Isobaric Interferences. Ionization suppression due to matrix effects was essentially eliminated by removing major cations (e.g., Ca^{2+} , Na^+) with the ion exchange column chromatography. However, isotopes of chromium and titanium as well as interfering polyatomic ions ClO^+ , ArC^+ , and PO^+ make the determination of V by isotope dilution especially complex (Table 4).

Elimination of ClO^+ Interference. The ClO^+ ion at m/z 51 is the primary interference in vanadium-51 measurement by ICPMS and prevents accurate determination in matrices with Cl levels as low as 0.01%. Sample introduction with ETV effectively eliminates the ClO^+ isobaric interference in solutions of 0.01–0.5 N HCl (Table 5). By preheating the sample in a graphite furnace before addition to the plasma, the volatile Cl is removed

⁽⁴³⁾ Heumann, K. G. Isotope Dilution Mass Spectrometry. In *Inorganic Mass Spectrometry*, Adams, F., Gijbels, R., Grieken, R., Eds.; J. Wiley and Sons: New York, 1988; pp 301–376.

New York, 1988; pp 301–376. (44) Adriaens, A. G.; Kelly, W. R.; Adams, R. C. *Anal. Chem.* **1993**, *65*, 660–

⁽⁴⁵⁾ Russ, G. P.; Bazan, J. M. Spectrochim. Acta 1987, 42B, 49-62.

⁽⁴⁶⁾ Gregoire, D. C. Anal. Chem. 1987, 59, 2479-2484.

⁽⁴⁷⁾ Jarvis, K. E.; Gray, A. L.; Houk, R. S. Handbook of inductively coupled plasma mass spectrometry, Blackie and Son: Glasgow, 1992; p 380.

Table 5. Evaluation of the Importance of the CIO⁺ Interference at *m/z* 51 as a Function of HCI Concentration Using Standard Pneumatic Nebulization Compared with Electrothermal Vaporization To Introduce the Sample^a

ETV at $m/z 51^b$ (counts/s)		
)		
i		

^a Volatilization of Cl in the graphite furnace before analysis eliminates the background due to ClO⁺, even at high concentrations of Cl. ^b Blank corrected values.

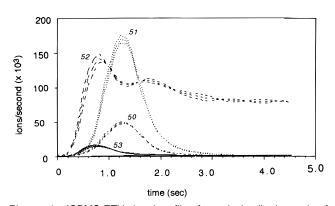


Figure 1. ICPMS-ETV signal profile of a typical spiked sample of foraminiferal calcite following cleaning and ion exchange. Shown are profiles of intensity at m/z 50, 51, 52, and 53. Peak intensity for m/z 49 is typically 1000 counts/s at 1.4 s and is at baseline in this figure. Note the high Cr-52 and 40 Ar 12 C peaks at 0.8 and 1.7 s, respectively, at m/z. Duplicate analyses and the averaged profile are shown.

and the ClO $^+$ interference eliminated. The pyrolysis temperature and time were optimized for complete removal of Cl with no loss of V. Based on these criteria, the optimum pyrolysis condition was found to be 1000 °C for 20 s. Even at only 200 °C for 5 s, 90% of the Cl was removed.

Titanium and Chromium-50 Correction. Since both Ti and Cr have isotopes at m/z 50 (relative abundances of 5.31 and 4.31%, respectively), they interfere with the ⁵⁰V measurement (Table 4). While efforts were made to separate these elements during the ion exchange process, measurable quantities of both were still present. ⁴⁹Ti (5.51%) was used to correct for the ⁵⁰Ti. Titanium isotopes at mass 48 and 46 were avoided due to the potential interference with calcium; ⁴⁷Ti has a higher and more variable background due to the PO+ ion. 50Ti accounts for 2-4% of the ⁵⁰V signal. The correction for the interference due to ⁵⁰Cr was more complex since isobaric elements exist at all the other chromium isotopes: ⁵²Cr (83.76%), ⁵³Cr (9.55%), and ⁵⁴Cr (2.38%). ⁵⁴Cr was not a viable alternative because of the high background from ⁴⁰Ar¹⁴N and the isobaric overlap with ⁵⁴Fe, which is not easily corrected. Polyatomic ions resulting from carbon argides at m/z52 (40 Ar 12 C) and 53 (40 Ar 13 C) overlap with 52 Cr and 53 Cr. They are difficult to correct for due to their very broad signal pulse, which is especially prominent at m/z 52 in Figure 1. Peak height measurements rather than integrated signals were used for the determination of Cr, which significantly reduces the background counts due to these carbon argides.³⁸ Measurement of ⁵³Cr was used to estimate ⁵⁰Cr because the intensity of the Ar¹³C peak is 2 orders of magnitude lower than the Ar¹²C peak (13 C = 0.011 12 C).

A value of 400-800 counts/s for Ar¹³C based on 1.1% of the Ar¹²C peak was subtracted from the total counts at m/z 53. The ⁵⁰Cr content was calculated from the resulting ⁵³Cr and used to correct the ⁵⁰V. Since Cr is vaporized off the graphite tube more readily than vanadium (0.8 s vs 1.3 s; Figure 1), and because peak height measurements are used, the Cr intensity is reduced to 60% of its maximum value at the vanadium peak. Thus, only 0.60 times the intensity of ⁵⁰Cr was used in the final correction. The equation used for the complete correction is given by

50
V = m/z 50 - 0.6(50 Cr) - (50 Ti)
= m/z 50 - 0.6(4.31/9.55)(m/z 53 - 40 Ar 13 C) - (5.34/5.51)(49 Ti)

 50 Cr levels typically ranged from 5 to 15% of the total signal at m/z 50. Figure 1 is a characteristic signal profile that illustrates the relative magnitude of vanadium isotope signals and interferences of a typical spiked sample. This sample contained an estimated 60 pg of natural vanadium and was spiked with 20 pg of 50 V and 25 pg of 51 V.

Detection Limit. The detection limit for vanadium by ETV-ICPMS without the isotope spike was 0.3 pg of V, based on the mass of vanadium that gives a response equivalent to 3 times the standard deviation of the blank at m/z 51. In this case, the blank was the response obtained from the vaporization of 33 μ L of 0.1 N HNO₃ and 20 μ L of the Pd modifier (n=15).

Using isotope dilution, the isobaric interferences at m/z 50 from the chromium and titanium result in a detection limit which is significantly higher. This was determined by treating the enriched isotope spike as a blank and measuring the V concentration on 0.1 mL aliquots. By definition, the vanadium concentration of the spike with no added sample, determined by isotope dilution, should be zero. The detection limit, defined as 3 times the standard deviation of the calculated value of the spike, was 6 pg of V. This higher value is due to the variability of the Cr and Ti levels, occasional ArC peaks from the graphite tube much higher than average, and the uncertainty in the correction for these isobaric interferences.

Sensitivity, Precision, and Accuracy of Sample Analysis. Sensitivity varied from run to run, depending on the condition of the cones, Ar flow, tuning of ion lenses, and other instrument conditions. Over a 1 year period, response at m/z 51 averaged 2500 counts s^{-1} pg^{-1} of V ($\pm 25\%$). An important matrix effect was observed that resulted in significant signal suppression above 1 mM Ca^{2+} . At 10 mM Ca, the vanadium intensity was 65% lower than that for a Ca-free matrix; at 100 mM Ca, signal suppression was 90%. The levels of Ca^{2+} in an untreated sample were ~ 0.5 M. Similar suppression was observed for a seawater matrix, necessitating the ion exchange to remove major cations, including Na^{+} and Ca^{2+} , from the sample. Following this separation step, Ca^{2+} levels in the dissolved calcite samples were < 10 nmol/L as determined by ICPMS at m/z 43.

Precision. Using 0.1 g of the dissolved foraminiferal standard, or the equivalent of about 2.5 mg of calcite after cleaning, the precision of the vanadium determination by ETV-ICPMS within one run was $\pm 3\%$ (1 σ ; n=10).

Accuracy. We assessed the accuracy of the analysis by two independent tests: standard addition and comparison to other methods for V determination. The first check was performed by

o i ·

Table 6. Intercomparison of Methods for Determination of V in Calcite

method	V (pg of V/g of CaCO ₃)	amount used (mg of CaCO ₃)	precision (%)
ID-ICPMS-ETV	370	2.5	$^{\pm 3}_{\pm 5}_{\pm 4}$
ID-TIMS	360	25	
APDC GFAAS	339	300	

adding 10-1000 pg of vanadium to $250~\mu L$ aliquots of the ultrapure CaCO₃ standard, which is equivalent to 7.5 mg of CaCO₃ or a typical foraminifera sample. These standards were processed according to the described method and analyzed by ID-ETV-ICPMS. Final values were corrected for background levels of V in the CaCO₃ standard, which was independently quantified by isotope dilution to be 14 pg of V/250 μL aliquot. A comparison of 10-1000 pg of V standard added to calcite versus vanadium of these samples, as determined by ID-ETV-ICPMS, resulted in a near perfect correlation. Vanadium added compared to the values as calculated by this method resulted in a linear correlation with a slope of 0.998 and $r^2=1.000$, which demonstrates that the method is accurate and linear over this range of values. Based on the instrument characteristics of the ICPMS, the method should be linear to levels at least 5 times greater, or 5000 pg.

Another indication of the accuracy of the measurement is to analyze the vanadium concentration of the standard foraminifera mixture by several different and independent methods: ID-ETV-ICPMS, ID-TIMS, and coprecipitation with APDC and graphite furnace atomic absorption spectrophotometry. As the results shown in Table 6 indicate, the two isotope dilution methods agree to within the precision of the measurement. The value for the

GFAAS determination was 10% lower than that obtained with the other methods, possibly due to incomplete recovery of V by the APDC.

CONCLUSIONS

Isotope dilution inductively coupled plasma mass spectrometry using electrothermal vaporization to remove isobaric interference at m/z 51 provides an extremely sensitive tool for the determination of V at picogram levels. This technique is sufficiently sensitive to allow the measurement of V in foraminiferal calcite, which has the potential to quantify changes in the redox state of the deep ocean over glacial-interglacial cycles. This method can be directly applied to the determination of V in other calcareous biogenic phases, including coralline aragonite and pteropods. Assuming that these living organisms incorporate vanadium into their skeleton in proportion to the vanadium concentration of seawater, foraminifera could also be useful paleoceanographic tools over different time scales and in different depositional environments. The technique is useful for the determination of V in seawater where sample volumes are limited, such as in benthic chamber experiments and porewater studies. The use of ⁵⁰V as an isotope tracer eliminates any inaccuracy due to variable recoveries during sample processing. Analysis of other environmental and biological samples characterized by ultralow levels of vanadium could also be facilitated by this method.

Received for review August 22, 1995. Accepted October 24, 1995.[®]

AC950859C

 $^{^{\}otimes}$ Abstract published in Advance ACS Abstracts, December 1, 1995.